

2023-12-21

URGENT FIELD SAFETY NOTICE

- Manufacturer SRN:** DE-MF-000020091
- FSCA Reference:** 946519 HLS – potentially compromised sterile barrier
- FSN Type:** New
- Affected Product:** Refer to Annex I
- Unique Device Identifier(s) (UDI-DI):** Refer to Annex I
- Affected Batch No.:** Refer to Annex I
- For Attention of:** Users of the medical devices listed in Annex I

Dear valued customer,

Maquet Cardiopulmonary GmbH (MCP)/Getinge is notifying you about an issue related to the HLS Sets potentially having a compromised sterile barrier of the Emergency Priming Line.

The HLS Set is intended for use in an extracorporeal circulation for cardiac support and/or pulmonary support. The Emergency Priming Line is for removing air in the blue line which has not yet reached the HLS module Set (Figure 1) and is therefore used in emergency situations only.

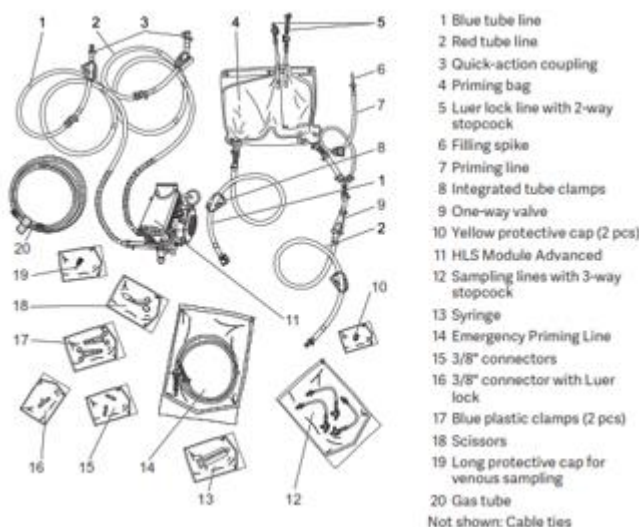


Figure 1: Blue tube line (1), Emergency Priming Line (14)

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

Description of the Problem

During retrospective internal testing of HLS Set packaging, MCP/Getinge discovered narrow channels in the seal of the sterile bags used to package the Emergency Priming Line (Figure 1 Item 14). MCP/Getinge receives this sterile bag from a supplier that performs tests to verify sterility. Incomplete sealing can potentially lead to a breach of the sterile barrier.

The Emergency Priming Line is first packaged in the above described sterile bag, which is itself packaged in the sterile Intellipack tray with the other HLS Set components. Thus, even with the potentially compromised sterile bag, the Emergency Priming Line is sterile when an HLS Set is first opened. Therefore, Emergency Priming Lines affected by this Field Safety Notice may only be used immediately after opening an HLS Set or else disposed of, if not needed. If an Emergency Priming Line is needed at a later time, a new HLS Set needs to be opened to obtain a sterile Emergency Priming Line.

To gauge the scope of affected products MCP/Getinge re-tested all lots of the supplier-provided packaging used for products that are still within shelf life. This testing demonstrated that three supplier batches of sterile bags were affected. Therefore, this Field Safety Notice is limited to products that contain sterile bags from the three affected supplier batches.

Hazardous situation / Risk to Health

- MCP/Getinge's Health Hazard Evaluation (HHE) determined that a breach of the sterile barrier of the Emergency Priming Line could expose patients to pathogenic agents.

This hazardous situation could result in the following potential harms (for further information, please refer to Annex II):

- Inflammation (medium risk)
- Infection (medium risk)
- Sepsis (medium risk)

There have been no adverse events or customer complaints reported related to the issue described above.

Corrective Action:

- **As described in the HLS Set Instructions for Use (IFU), device users should always have a backup HLS Set in stock and available during therapy.** Please check your inventory to ensure you have a backup HLS Set.
1. If you have a backup HLS Set and need to use the Emergency Priming Line:
 - Dispose of the Emergency Priming Line from the affected HLS Set.
 - Use the Emergency Priming Line from the backup HLS Set.
 - The remaining components of the backup HLS Set must be disposed of after extraction of the Emergency Priming Line.
 - Please reach out to us for the issuance of a replacement HLS Set or credit note.
 2. If you do not have a backup HLS Set:
 - Do not use the HLS Set for therapy. Please contact your local customer service to request a backup before using the product. We will send you a backup set within 5-7 days upon request.
 - For future needs, please contact us to discuss receiving an additional backup set.

Action to be taken by the user:

- | | |
|---|--|
| <input checked="" type="checkbox"/> Identify Device | <input type="checkbox"/> Quarantine Device |
| <input checked="" type="checkbox"/> Return Device | <input type="checkbox"/> Destroy Device |

Details of the further action(s):

- According to our post-market surveillance documentation, you may have products affected by this action. Please examine your inventory immediately to determine, if you have the affected HLS in your inventory.
- If a product is already in use, it should remain in use due to increased potential risk when disconnecting the product during ongoing therapy.
- If you need to use an Emergency Priming Line, you must dispose of the Emergency Priming Line from any affected products and take a new Emergency Priming Line from a backup unopened HLS Set.
- Please be aware that the components supplied with an HLS Set must be consumed **immediately after opening the Tyvek lid** of the HLS Set tray. **The remaining components of a newly opened HLS set must be disposed of** after extraction of the Emergency Priming Line.
- Please **always** report any adverse events, e.g., infections potentially related to the affected products, to your Getinge representative.
- Whether or not you have affected product(s), duly fill out the enclosed Customer Response Form / Letter of Acknowledgement and return it to your local Getinge representative by **January 17, 2024**, the latest. Please give **FSCA-946519** as reference in the subject line of your email. Upon return, Getinge will initiate the necessary measures according to the two options above.

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

Action to be taken by the manufacturer:

- Product Removal
- Software upgrade
- Other
- On-site device modification/ inspection
- IFU or labelling change
- None

- Inform all customers possessing the affected products **promptly** about this Field Action by sending the Field Safety Notice.
- Depending on the feedback from the customer:
 - If a customer uses the Emergency Priming Line from a backup HLS Set:
 - Provide the customer with a replacement HLS Set or credit note
 - If a customer has one leftover affected product with no backup product:
 - Provide the customer with a replacement HLS Set or arrange the return of the affected product
- Getinge is implementing additional corrective and preventive measures with the supplier in order to prevent the issue from happening again.
- Verification testing for integrity of the sterile seal is now performed by MCP/Getinge in-house before the material is used in production.

Enclosed documents:

- Customer response form
- Annex I List of affected products
- Annex II Further information regarding Hazardous situation, Harms and Risk Levels

Transmission of the Field Safety Notice

- Please ensure in your organization that all users of the above-mentioned products and any other persons that need to be informed are made aware of this Urgent Field Safety Notice.
- Please transfer this notice to other organizations that may be impacted by this action.
- If you have given the affected product(s) to third parties, please forward a copy of this information or inform the MCP/Getinge contact person indicated below.

Please remain aware of this notice and its accompanying actions for an appropriate period to ensure effectiveness of the corrective action

We sincerely apologize for any inconvenience this may cause you and will do our utmost to carry through this action as swiftly as possible.

As required, we have provided this notification to the necessary Regulatory Agencies.

Should you have questions or require additional information, please contact your local Getinge representative, or send an e-mail to FSCA.cp@getinge.com.

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

Sincerely,

Managing Director

Signature: *Dieter Engel*

Electronically signed by: Dieter Engel
Reason: I approve this document.
Date: Dec 21, 2023 13:24 GMT+1

Email: dieter.engel@getinge.com

**Person Responsible for Regulatory
Compliance (PRRC)
(on behalf of the PRRC)**

Signature: *Alexander Bernhardt*

Electronically signed by: Alexander
Bernhardt
Reason: I approve this document.
Date: Dec 21, 2023 12:56 GMT+1

Email: alexander.bernhardt@getinge.com

Contact details of manufacturer

Tom Peters
Maquet Cardiopulmonary GmbH
Kehler Str. 31
76437 Rastatt
GERMANY
Phone: +49 7222 932 - 0
Email: FSCA.cp@getinge.com

CUSTOMER RESPONSE FORM

FSCA Reference: 946519 HLS– potentially compromised sterile barrier

Affected Product: Refer to Annex I

Affected Batch No.: Refer to Annex I

Please send this form at the latest by **January 24, 2024**, to your local Getinge representative.

By completing this document and signing it, I acknowledge that I have read and understand the following associated points:

- I have read and understand this Field Safety Notice for affected product HLS. We will take action as soon as possible according to given instructions.
- I confirm that I have distributed this Field Safety Notice to the affected personal.

I do not have any HLS in my inventory.

I have following HLS in my inventory and decided for the following option:

| Article No. | Description | Batch No. | Quantity |
|-------------|-------------|-----------|----------|
| | | | |
| | | | |
| | | | |
| | | | |

Your Comments:

Country

Hospital / Clinic (full address)

Date

Name (Function)

Signature

Please return the completed form to your local Getinge representative by email [enter local Getinge mail address](#) or via post [enter local Getinge address](#) or FAX.

Annex I List of affected products

This Annex I List of affected products is considered a supplementary attachment to the 946519 Field Safety Notice.

| COUNTRY | ITEM | UDI | LOT OR SERIAL # | ITEM DESCRIPTION |
|-----------|-----------|---------------|-----------------|---|
| Argentina | 701069073 | 4058863005744 | 3000281237 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Argentina | 701069073 | 4058863005744 | 3000327125 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Australia | 701069073 | 4058863005744 | 3000281883 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Australia | 701069073 | 4058863005744 | 3000283435 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Australia | 701069065 | 4058863300238 | 3000301224 | BEQ-HLS 7050-CA #SHLS Set Advanced 7.0 |
| Australia | 701069065 | 4058863300238 | 3000317530 | BEQ-HLS 7050-CA #SHLS Set Advanced 7.0 |
| Australia | 701069073 | 4058863005744 | 3000317548 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Australia | 701069073 | 4058863005744 | 3000319746 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Australia | 701069065 | 4058863300238 | 3000324493 | BEQ-HLS 7050-CA #SHLS Set Advanced 7.0 |
| Australia | 701069073 | 4058863005744 | 3000319774 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Australia | 701069073 | 4058863005744 | 3000323079 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Australia | 701069076 | 4058863078298 | 3000323402 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Australia | 701069073 | 4058863005744 | 3000324505 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Australia | 701069073 | 4058863005744 | 3000325553 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Australia | 701069077 | 4058863076355 | 3000327602 | BEQ-HLS 5050 USA #SHLS Set Advanced 5.0 |
| Australia | 701069078 | 4058863080383 | 3000323082 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| Australia | 701069078 | 4058863080383 | 3000323084 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| Australia | 701069078 | 4058863080383 | 3000324496 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| Australia | 701069078 | 4058863080383 | 3000324499 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| Australia | 701069078 | 4058863080383 | 3000327120 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|-----------|-----------|---------------|------------|---|
| Australia | 701069078 | 4058863080383 | 3000328404 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| Austria | 701069073 | 4058863005744 | 3000303593 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Austria | 701069073 | 4058863005744 | 3000316207 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Austria | 701069073 | 4058863005744 | 3000303595 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Austria | 701069073 | 4058863005744 | 3000298132 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Austria | 701069073 | 4058863005744 | 3000318848 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Austria | 701069073 | 4058863005744 | 3000319744 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Austria | 701069073 | 4058863005744 | 3000319746 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Austria | 701069073 | 4058863005744 | 3000319774 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Austria | 701069073 | 4058863005744 | 3000323079 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Austria | 701069073 | 4058863005744 | 3000325552 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Austria | 701069073 | 4058863005744 | 3000325553 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Austria | 701069073 | 4058863005744 | 3000324503 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Austria | 701069073 | 4058863005744 | 3000327125 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Austria | 701069073 | 4058863005744 | 3000318847 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Austria | 701069073 | 4058863005744 | 3000336493 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Bahrain | 701069073 | 4058863005744 | 3000282262 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Bahrain | 701069073 | 4058863005744 | 3000316207 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Bahrain | 701069073 | 4058863005744 | 3000324505 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Belarus | 701069073 | 4058863005744 | 3000317547 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Belarus | 701069073 | 4058863005744 | 3000317548 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Belarus | 701069073 | 4058863005744 | 3000323079 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Belarus | 701069073 | 4058863005744 | 3000324505 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|---------|-----------|---------------|------------|---|
| Belgium | 701069073 | 4058863005744 | 3000301225 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Belgium | 701069073 | 4058863005744 | 3000303593 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Belgium | 701069073 | 4058863005744 | 3000316207 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Belgium | 701069073 | 4058863005744 | 3000303595 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Belgium | 701069073 | 4058863005744 | 3000317548 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Belgium | 701069073 | 4058863005744 | 3000318848 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Belgium | 701069073 | 4058863005744 | 3000319744 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Belgium | 701069073 | 4058863005744 | 3000319746 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Belgium | 701069073 | 4058863005744 | 3000319774 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Belgium | 701069073 | 4058863005744 | 3000323079 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Belgium | 701069073 | 4058863005744 | 3000324505 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Belgium | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Belgium | 701069073 | 4058863005744 | 3000327125 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Belgium | 701069073 | 4058863005744 | 3000318847 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Belgium | 701069073 | 4058863005744 | 3000336493 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Canada | 701069065 | 4058863300238 | 3000301223 | BEQ-HLS 7050-CA #SHLS Set Advanced 7.0 |
| Canada | 701069065 | 4058863300238 | 3000301224 | BEQ-HLS 7050-CA #SHLS Set Advanced 7.0 |
| Canada | 701069065 | 4058863300238 | 3000317530 | BEQ-HLS 7050-CA #SHLS Set Advanced 7.0 |
| Canada | 701069063 | 4058863304632 | 3000313299 | BO-HLS 7050-CA #S/ HIT Set Advanced 7.0 |
| Canada | 701069063 | 4058863304632 | 3000319771 | BO-HLS 7050-CA #S/ HIT Set Advanced 7.0 |
| Canada | 701069068 | 4058863304625 | 3000323405 | BEQ-HLS 5050-CA #SHLS Set Advanced 5.0 |
| Canada | 701069065 | 4058863300238 | 3000324493 | BEQ-HLS 7050-CA #SHLS Set Advanced 7.0 |
| Canada | 701069068 | 4058863304625 | 3000329236 | BEQ-HLS 5050-CA #SHLS Set Advanced 5.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|------------|-----------|---------------|------------|--|
| Canada | 701069065 | 4058863300238 | 3000337837 | BEQ-HLS 7050-CA #SHLS Set Advanced 7.0 |
| Chile | 701069073 | 4058863005744 | 3000316210 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Chile | 701069073 | 4058863005744 | 3000319744 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Chile | 701069073 | 4058863005744 | 3000325551 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Chile | 701069073 | 4058863005744 | 3000327125 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| China | 701069073 | 4058863005744 | 3000303593 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| China | 701069073 | 4058863005744 | 3000316207 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| China | 701069076 | 4058863078298 | 3000311768 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| China | 701069073 | 4058863005744 | 3000317547 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| China | 701069073 | 4058863005744 | 3000318848 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| China | 701069073 | 4058863005744 | 3000319744 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| China | 701069073 | 4058863005744 | 3000319746 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| China | 701069073 | 4058863005744 | 3000319774 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| China | 701069073 | 4058863005744 | 3000323079 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| China | 701069076 | 4058863078298 | 3000323403 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| China | 701069073 | 4058863005744 | 3000325551 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| China | 701069076 | 4058863078298 | 3000327124 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| China | 701069073 | 4058863005744 | 3000325553 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| China | 701069073 | 4058863005744 | 3000325554 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| China | 701069073 | 4058863005744 | 3000318847 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| China | 701069073 | 4058863005744 | 3000336494 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Colombia | 701069073 | 4058863005744 | 3000316210 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Costa Rica | 701069073 | 4058863005744 | 3000301226 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|----------------|-----------|---------------|------------|------------------------------------|
| Costa Rica | 701069073 | 4058863005744 | 3000316207 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Costa Rica | 701069073 | 4058863005744 | 3000316210 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Costa Rica | 701069073 | 4058863005744 | 3000319774 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Costa Rica | 701069073 | 4058863005744 | 3000325554 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Croatia | 701069073 | 4058863005744 | 3000301227 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Croatia | 701069073 | 4058863005744 | 3000298132 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Croatia | 701069073 | 4058863005744 | 3000324503 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Cyprus | 701069073 | 4058863005744 | 3000316210 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Cyprus | 701069073 | 4058863005744 | 3000319744 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Cyprus | 701069073 | 4058863005744 | 3000319774 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Cyprus | 701069073 | 4058863005744 | 3000324504 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Cyprus | 701069073 | 4058863005744 | 3000325554 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Cyprus | 701069073 | 4058863005744 | 3000325551 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Czech Republic | 701069073 | 4058863005744 | 3000282259 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Czech Republic | 701069073 | 4058863005744 | 3000282260 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Czech Republic | 701069073 | 4058863005744 | 3000303593 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Czech Republic | 701069073 | 4058863005744 | 3000316207 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Czech Republic | 701069073 | 4058863005744 | 3000303595 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Czech Republic | 701069073 | 4058863005744 | 3000317550 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Czech Republic | 701069073 | 4058863005744 | 3000318848 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Czech Republic | 701069073 | 4058863005744 | 3000319744 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Czech Republic | 701069073 | 4058863005744 | 3000319746 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Czech Republic | 701069073 | 4058863005744 | 3000319774 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|----------------|-----------|---------------|------------|------------------------------------|
| Czech Republic | 701069073 | 4058863005744 | 3000323079 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Czech Republic | 701069073 | 4058863005744 | 3000325551 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Czech Republic | 701069073 | 4058863005744 | 3000325554 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Czech Republic | 701069073 | 4058863005744 | 3000324503 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Czech Republic | 701069073 | 4058863005744 | 3000324504 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Czech Republic | 701069073 | 4058863005744 | 3000327125 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Czech Republic | 701069073 | 4058863005744 | 3000318847 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Czech Republic | 701069073 | 4058863005744 | 3000336494 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Denmark | 701069073 | 4058863005744 | 3000281237 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Denmark | 701069073 | 4058863005744 | 3000281885 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Denmark | 701069073 | 4058863005744 | 3000303593 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Denmark | 701069073 | 4058863005744 | 3000316207 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Denmark | 701069073 | 4058863005744 | 3000303595 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Denmark | 701069073 | 4058863005744 | 3000298132 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Denmark | 701069073 | 4058863005744 | 3000318848 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Denmark | 701069073 | 4058863005744 | 3000319774 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Denmark | 701069073 | 4058863005744 | 3000323079 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Denmark | 701069076 | 4058863078298 | 3000327124 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Denmark | 701069073 | 4058863005744 | 3000317551 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Denmark | 701069073 | 4058863005744 | 3000336494 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Egypt | 701069073 | 4058863005744 | 3000282259 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Egypt | 701069073 | 4058863005744 | 3000282194 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Finland | 701069073 | 4058863005744 | 3000282259 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|---------|-----------|---------------|------------|--------------------------------------|
| Finland | 701069073 | 4058863005744 | 3000303594 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Finland | 701069073 | 4058863005744 | 3000314516 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Finland | 701069073 | 4058863005744 | 3000317533 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Finland | 701069073 | 4058863005744 | 3000303595 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Finland | 701069073 | 4058863005744 | 3000298132 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Finland | 701069073 | 4058863005744 | 3000318849 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Finland | 701069073 | 4058863005744 | 3000319745 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Finland | 701069073 | 4058863005744 | 3000319747 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Finland | 701069073 | 4058863005744 | 3000319774 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Finland | 701069073 | 4058863005744 | 3000323079 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Finland | 701069073 | 4058863005744 | 3000324503 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Finland | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Finland | 701069073 | 4058863005744 | 3000327125 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Finland | 701069073 | 4058863005744 | 3000336494 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069073 | 4058863005744 | 3000281885 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069073 | 4058863005744 | 3000282194 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069073 | 4058863005744 | 3000303594 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069073 | 4058863005744 | 3000314516 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069073 | 4058863005744 | 3000317533 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069073 | 4058863005744 | 3000303595 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069073 | 4058863005744 | 3000317549 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069073 | 4058863005744 | 3000318845 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069083 | 4058863020082 | 3000311765 | BO-HLS 7050 #S/ HIT Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|------------------|-----------|---------------|------------|------------------------------------|
| France | 701069073 | 4058863005744 | 3000318849 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069073 | 4058863005744 | 3000319745 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069073 | 4058863005744 | 3000317548 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069073 | 4058863005744 | 3000319747 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069073 | 4058863005744 | 3000319775 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069073 | 4058863005744 | 3000324502 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069076 | 4058863078298 | 3000323402 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| France | 701069073 | 4058863005744 | 3000324503 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069076 | 4058863078298 | 3000327124 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| France | 701069073 | 4058863005744 | 3000325553 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069073 | 4058863005744 | 3000325552 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069073 | 4058863005744 | 3000325554 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069076 | 4058863078298 | 3000323403 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| France | 701069073 | 4058863005744 | 3000327125 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069073 | 4058863005744 | 3000318847 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| France | 701069076 | 4058863078298 | 3000327600 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| France | 701069073 | 4058863005744 | 3000336493 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| French Polynesia | 701069073 | 4058863005744 | 3000318845 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| French Polynesia | 701069073 | 4058863005744 | 3000319745 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000281237 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000282260 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000283435 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|---------|-----------|---------------|------------|--------------------------------------|
| Germany | 701069073 | 4058863005744 | 3000301227 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000301226 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000303593 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000303594 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000316209 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069076 | 4058863078298 | 3000311768 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Germany | 701069073 | 4058863005744 | 3000303595 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000317550 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000317551 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000318848 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000317548 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000319744 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000319746 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000319745 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000319775 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069083 | 4058863020082 | 3000311765 | BO-HLS 7050 #S/ HIT Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000324501 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069076 | 4058863078298 | 3000323404 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Germany | 701069073 | 4058863005744 | 3000325552 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000324504 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069076 | 4058863078298 | 3000327124 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Germany | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000325556 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|-----------|-----------|---------------|------------|------------------------------------|
| Germany | 701069076 | 4058863078298 | 3000327600 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Germany | 701069073 | 4058863005744 | 3000301225 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000327125 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000318847 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000336493 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000316210 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000324502 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000324505 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000325551 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000325553 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069076 | 4058863078298 | 3000323403 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Germany | 701069073 | 4058863005744 | 3000336494 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Germany | 701069073 | 4058863005744 | 3000336495 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Greece | 701069073 | 4058863005744 | 3000301227 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Greece | 701069073 | 4058863005744 | 3000317533 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Greece | 701069073 | 4058863005744 | 3000298132 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Greece | 701069073 | 4058863005744 | 3000319744 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hong Kong | 701069073 | 4058863005744 | 3000283435 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hong Kong | 701069073 | 4058863005744 | 3000282260 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hong Kong | 701069073 | 4058863005744 | 3000301226 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hong Kong | 701069073 | 4058863005744 | 3000303594 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hong Kong | 701069073 | 4058863005744 | 3000317533 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hong Kong | 701069073 | 4058863005744 | 3000316210 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|-----------|-----------|---------------|------------|------------------------------------|
| Hong Kong | 701069073 | 4058863005744 | 3000318849 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hong Kong | 701069073 | 4058863005744 | 3000319745 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hong Kong | 701069073 | 4058863005744 | 3000319747 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hong Kong | 701069073 | 4058863005744 | 3000319775 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hong Kong | 701069073 | 4058863005744 | 3000324501 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hong Kong | 701069073 | 4058863005744 | 3000324503 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hong Kong | 701069073 | 4058863005744 | 3000325554 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hong Kong | 701069073 | 4058863005744 | 3000318847 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hungary | 701069073 | 4058863005744 | 3000301227 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hungary | 701069073 | 4058863005744 | 3000317533 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hungary | 701069073 | 4058863005744 | 3000317549 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hungary | 701069073 | 4058863005744 | 3000298132 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hungary | 701069073 | 4058863005744 | 3000318849 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hungary | 701069073 | 4058863005744 | 3000319745 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hungary | 701069073 | 4058863005744 | 3000319747 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hungary | 701069076 | 4058863078298 | 3000323402 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Hungary | 701069073 | 4058863005744 | 3000324503 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hungary | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hungary | 701069073 | 4058863005744 | 3000319746 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Hungary | 701069073 | 4058863005744 | 3000336493 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| India | 701069073 | 4058863005744 | 3000282262 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| India | 701069073 | 4058863005744 | 3000317533 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| India | 701069073 | 4058863005744 | 3000316210 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|-----------|-----------|---------------|------------|------------------------------------|
| India | 701069073 | 4058863005744 | 3000324503 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| India | 701069073 | 4058863005744 | 3000325553 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| India | 701069073 | 4058863005744 | 3000318847 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Indonesia | 701069073 | 4058863005744 | 3000319747 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Iraq | 701069073 | 4058863005744 | 3000282261 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Iraq | 701069073 | 4058863005744 | 3000282262 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Ireland | 701069073 | 4058863005744 | 3000281237 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Ireland | 701069073 | 4058863005744 | 3000281883 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Ireland | 701069073 | 4058863005744 | 3000301226 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Ireland | 701069073 | 4058863005744 | 3000319745 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Ireland | 701069073 | 4058863005744 | 3000314516 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Israel | 701069073 | 4058863005744 | 3000303594 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Israel | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Italy | 701069073 | 4058863005744 | 3000281883 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Italy | 701069073 | 4058863005744 | 3000281885 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Italy | 701069073 | 4058863005744 | 3000303594 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Italy | 701069073 | 4058863005744 | 3000314516 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Italy | 701069073 | 4058863005744 | 3000316209 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Italy | 701069076 | 4058863078298 | 3000311768 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Italy | 701069073 | 4058863005744 | 3000317549 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Italy | 701069073 | 4058863005744 | 3000298132 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Italy | 701069073 | 4058863005744 | 3000318849 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Italy | 701069073 | 4058863005744 | 3000319745 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|----------------------------------|-----------|---------------|------------|--|
| Italy | 701069073 | 4058863005744 | 3000319747 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Italy | 701069073 | 4058863005744 | 3000319776 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Italy | 701069073 | 4058863005744 | 3000324502 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Italy | 701069076 | 4058863078298 | 3000323403 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Italy | 701069073 | 4058863005744 | 3000324503 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Italy | 701069076 | 4058863078298 | 3000327124 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Italy | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Italy | 701069073 | 4058863005744 | 3000317533 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Italy | 701069073 | 4058863005744 | 3000325556 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Italy | 701069073 | 4058863005744 | 3000324501 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Italy | 701069073 | 4058863005744 | 3000327125 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Italy | 701069073 | 4058863005744 | 3000318847 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Italy | 701069083 | 4058863020082 | 3000311765 | BO-HLS 7050 #S/ HIT Set Advanced 7.0 |
| Italy | 701069073 | 4058863005744 | 3000336494 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Japan | 701077943 | 4058863308777 | 3000324291 | BE-HLS 7050#Non-CE SHLS Set Advanced 7.0 |
| Japan | 701077943 | 4058863308777 | 3000324290 | BE-HLS 7050#Non-CE SHLS Set Advanced 7.0 |
| Jordan | 701069073 | 4058863005744 | 3000282262 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Kazakhstan | 701069073 | 4058863005744 | 3000317547 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Korea, Republic of (South Korea) | 701069073 | 4058863005744 | 3000301225 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Korea, Republic of (South Korea) | 701069073 | 4058863005744 | 3000303594 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Korea, Republic of (South Korea) | 701069073 | 4058863005744 | 3000317533 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|----------------------------------|-----------|---------------|------------|------------------------------------|
| Korea, Republic of (South Korea) | 701069073 | 4058863005744 | 3000316210 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Korea, Republic of (South Korea) | 701069073 | 4058863005744 | 3000318849 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Korea, Republic of (South Korea) | 701069073 | 4058863005744 | 3000319745 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Korea, Republic of (South Korea) | 701069073 | 4058863005744 | 3000319747 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Korea, Republic of (South Korea) | 701069073 | 4058863005744 | 3000324501 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Korea, Republic of (South Korea) | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Kuwait | 701069073 | 4058863005744 | 3000282194 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Kuwait | 701069073 | 4058863005744 | 3000281237 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Kuwait | 701069073 | 4058863005744 | 3000303594 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Kuwait | 701069073 | 4058863005744 | 3000336494 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Luxembourg | 701069073 | 4058863005744 | 3000317533 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Luxembourg | 701069073 | 4058863005744 | 3000316210 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Luxembourg | 701069073 | 4058863005744 | 3000319776 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Malaysia | 701069073 | 4058863005744 | 3000301226 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Martinique | 701069073 | 4058863005744 | 3000314516 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Martinique | 701069073 | 4058863005744 | 3000318849 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Martinique | 701069073 | 4058863005744 | 3000319747 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Martinique | 701069073 | 4058863005744 | 3000324502 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Martinique | 701069076 | 4058863078298 | 3000323402 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Martinique | 701069073 | 4058863005744 | 3000324503 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|-------------|-----------|---------------|------------|------------------------------------|
| Martinique | 701069076 | 4058863078298 | 3000327124 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Martinique | 701069073 | 4058863005744 | 3000327125 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Martinique | 701069073 | 4058863005744 | 3000336493 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Mayotte | 701069073 | 4058863005744 | 3000319775 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Mayotte | 701069076 | 4058863078298 | 3000323402 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Mexico | 701069073 | 4058863005744 | 3000282261 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Mexico | 701069073 | 4058863005744 | 3000317533 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Mexico | 701069073 | 4058863005744 | 3000316210 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Mexico | 701069073 | 4058863005744 | 3000324502 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Mexico | 701069073 | 4058863005744 | 3000319747 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Mexico | 701069076 | 4058863078298 | 3000327124 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Netherlands | 701069073 | 4058863005744 | 3000301225 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Netherlands | 701069073 | 4058863005744 | 3000303594 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Netherlands | 701069076 | 4058863078298 | 3000311768 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Netherlands | 701069073 | 4058863005744 | 3000316209 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Netherlands | 701069073 | 4058863005744 | 3000317549 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Netherlands | 701069073 | 4058863005744 | 3000298132 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Netherlands | 701069073 | 4058863005744 | 3000318849 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Netherlands | 701069073 | 4058863005744 | 3000319745 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Netherlands | 701069073 | 4058863005744 | 3000319747 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Netherlands | 701069073 | 4058863005744 | 3000319775 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Netherlands | 701069073 | 4058863005744 | 3000324502 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Netherlands | 701069076 | 4058863078298 | 3000323402 | BE-HLS 5050 #SHLS Set Advanced 5.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|---------------|-----------|---------------|------------|--------------------------------------|
| Netherlands | 701069073 | 4058863005744 | 3000323079 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Netherlands | 701069076 | 4058863078298 | 3000323404 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Netherlands | 701069073 | 4058863005744 | 3000324504 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Netherlands | 701069076 | 4058863078298 | 3000327124 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Netherlands | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Netherlands | 701069073 | 4058863005744 | 3000316207 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Netherlands | 701069073 | 4058863005744 | 3000327125 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Netherlands | 701069073 | 4058863005744 | 3000318847 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Netherlands | 701069073 | 4058863005744 | 3000336494 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| New Caledonia | 701069073 | 4058863005744 | 3000317548 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| New Caledonia | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| New Caledonia | 701069073 | 4058863005744 | 3000318847 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| New Zealand | 701069073 | 4058863005744 | 3000324504 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Norway | 701069073 | 4058863005744 | 3000281883 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Norway | 701069073 | 4058863005744 | 3000303594 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Norway | 701069083 | 4058863020082 | 3000311765 | BO-HLS 7050 #S/ HIT Set Advanced 7.0 |
| Norway | 701069073 | 4058863005744 | 3000316209 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Norway | 701069073 | 4058863005744 | 3000317549 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Norway | 701069073 | 4058863005744 | 3000298132 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Norway | 701069073 | 4058863005744 | 3000318849 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Norway | 701069073 | 4058863005744 | 3000319745 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Norway | 701069073 | 4058863005744 | 3000319747 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Norway | 701069073 | 4058863005744 | 3000319776 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|-------------|-----------|---------------|------------|------------------------------------|
| Norway | 701069073 | 4058863005744 | 3000323079 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Norway | 701069076 | 4058863078298 | 3000323402 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Norway | 701069073 | 4058863005744 | 3000324503 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Norway | 701069073 | 4058863005744 | 3000325554 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Norway | 701069073 | 4058863005744 | 3000327125 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Norway | 701069073 | 4058863005744 | 3000318847 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Norway | 701069073 | 4058863005744 | 3000336494 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Oman | 701069073 | 4058863005744 | 3000282262 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Oman | 701069073 | 4058863005744 | 3000283435 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Oman | 701069073 | 4058863005744 | 3000324504 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Philippines | 701069073 | 4058863005744 | 3000303594 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Philippines | 701069073 | 4058863005744 | 3000317547 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Philippines | 701069073 | 4058863005744 | 3000319776 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Poland | 701069073 | 4058863005744 | 3000282194 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Poland | 701069073 | 4058863005744 | 3000281885 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Poland | 701069073 | 4058863005744 | 3000282258 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Poland | 701069073 | 4058863005744 | 3000282259 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Poland | 701069073 | 4058863005744 | 3000317533 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Poland | 701069073 | 4058863005744 | 3000317549 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Poland | 701069073 | 4058863005744 | 3000298132 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Poland | 701069073 | 4058863005744 | 3000318849 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Poland | 701069073 | 4058863005744 | 3000319745 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Poland | 701069073 | 4058863005744 | 3000319747 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|----------|-----------|---------------|------------|------------------------------------|
| Poland | 701069073 | 4058863005744 | 3000319776 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Poland | 701069073 | 4058863005744 | 3000323079 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Poland | 701069073 | 4058863005744 | 3000324504 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Poland | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Poland | 701069073 | 4058863005744 | 3000316207 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Poland | 701069073 | 4058863005744 | 3000318847 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Poland | 701069073 | 4058863005744 | 3000336494 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Portugal | 701069073 | 4058863005744 | 3000281237 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Portugal | 701069073 | 4058863005744 | 3000281883 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Portugal | 701069073 | 4058863005744 | 3000301226 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Portugal | 701069073 | 4058863005744 | 3000303594 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Portugal | 701069073 | 4058863005744 | 3000317533 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Portugal | 701069076 | 4058863078298 | 3000311768 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Portugal | 701069073 | 4058863005744 | 3000317549 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Portugal | 701069073 | 4058863005744 | 3000317551 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Portugal | 701069073 | 4058863005744 | 3000317550 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Portugal | 701069073 | 4058863005744 | 3000318849 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Portugal | 701069073 | 4058863005744 | 3000319745 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Portugal | 701069073 | 4058863005744 | 3000319747 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Portugal | 701069073 | 4058863005744 | 3000319776 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Portugal | 701069073 | 4058863005744 | 3000324502 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Portugal | 701069076 | 4058863078298 | 3000323402 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Portugal | 701069073 | 4058863005744 | 3000324504 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|----------|-----------|---------------|------------|------------------------------------|
| Portugal | 701069076 | 4058863078298 | 3000327124 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Portugal | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Portugal | 701069076 | 4058863078298 | 3000327600 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Portugal | 701069073 | 4058863005744 | 3000319744 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Portugal | 701069073 | 4058863005744 | 3000325552 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Portugal | 701069073 | 4058863005744 | 3000325551 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Portugal | 701069076 | 4058863078298 | 3000323404 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Portugal | 701069073 | 4058863005744 | 3000327125 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Portugal | 701069073 | 4058863005744 | 3000318847 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Portugal | 701069073 | 4058863005744 | 3000336493 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Qatar | 701069073 | 4058863005744 | 3000282194 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Qatar | 701069073 | 4058863005744 | 3000303594 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Qatar | 701069073 | 4058863005744 | 3000316207 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Qatar | 701069073 | 4058863005744 | 3000319776 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Réunion | 701069073 | 4058863005744 | 3000314516 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Réunion | 701069073 | 4058863005744 | 3000317533 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Réunion | 701069073 | 4058863005744 | 3000318845 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Réunion | 701069076 | 4058863078298 | 3000323402 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Réunion | 701069073 | 4058863005744 | 3000324503 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Réunion | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Réunion | 701069076 | 4058863078298 | 3000327124 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Réunion | 701069076 | 4058863078298 | 3000323403 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Romania | 701069073 | 4058863005744 | 3000301227 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|--------------|-----------|---------------|------------|------------------------------------|
| Romania | 701069073 | 4058863005744 | 3000317549 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Romania | 701069073 | 4058863005744 | 3000298132 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Romania | 701069073 | 4058863005744 | 3000318849 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Romania | 701069073 | 4058863005744 | 3000319745 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Romania | 701069073 | 4058863005744 | 3000319747 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Romania | 701069073 | 4058863005744 | 3000324504 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Romania | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Romania | 701069073 | 4058863005744 | 3000314516 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Romania | 701069076 | 4058863078298 | 3000323402 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Romania | 701069073 | 4058863005744 | 3000318847 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Romania | 701069073 | 4058863005744 | 3000336494 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Saudi Arabia | 701069073 | 4058863005744 | 3000282258 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Saudi Arabia | 701069073 | 4058863005744 | 3000283435 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Saudi Arabia | 701069073 | 4058863005744 | 3000282262 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Serbia | 701069073 | 4058863005744 | 3000317549 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Serbia | 701069073 | 4058863005744 | 3000298132 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Serbia | 701069073 | 4058863005744 | 3000318849 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Serbia | 701069073 | 4058863005744 | 3000319745 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Serbia | 701069073 | 4058863005744 | 3000319747 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Serbia | 701069073 | 4058863005744 | 3000324501 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Serbia | 701069073 | 4058863005744 | 3000324504 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Serbia | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Serbia | 701069073 | 4058863005744 | 3000316209 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|-----------|-----------|---------------|------------|------------------------------------|
| Singapore | 701069073 | 4058863005744 | 3000301226 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Singapore | 701069073 | 4058863005744 | 3000316210 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Singapore | 701069073 | 4058863005744 | 3000317547 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Singapore | 701069073 | 4058863005744 | 3000319774 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Singapore | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Slovakia | 701069073 | 4058863005744 | 3000301227 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Slovakia | 701069073 | 4058863005744 | 3000317549 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Slovakia | 701069073 | 4058863005744 | 3000316210 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Slovakia | 701069073 | 4058863005744 | 3000318849 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Slovakia | 701069073 | 4058863005744 | 3000319745 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Slovakia | 701069073 | 4058863005744 | 3000319747 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Slovakia | 701069073 | 4058863005744 | 3000319776 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Slovakia | 701069073 | 4058863005744 | 3000324502 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Slovakia | 701069073 | 4058863005744 | 3000324505 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Slovakia | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Slovakia | 701069073 | 4058863005744 | 3000323079 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Slovenia | 701069073 | 4058863005744 | 3000301227 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Slovenia | 701069073 | 4058863005744 | 3000317533 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Slovenia | 701069073 | 4058863005744 | 3000317549 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Slovenia | 701069073 | 4058863005744 | 3000318845 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Slovenia | 701069073 | 4058863005744 | 3000318849 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Slovenia | 701069073 | 4058863005744 | 3000319747 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Slovenia | 701069073 | 4058863005744 | 3000324502 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|--------------|-----------|---------------|------------|------------------------------------|
| Slovenia | 701069073 | 4058863005744 | 3000324504 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Slovenia | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Slovenia | 701069073 | 4058863005744 | 3000323079 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| South Africa | 701069073 | 4058863005744 | 3000303594 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069073 | 4058863005744 | 3000301225 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069073 | 4058863005744 | 3000303594 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069073 | 4058863005744 | 3000303593 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069073 | 4058863005744 | 3000314516 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069076 | 4058863078298 | 3000311768 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Spain | 701069073 | 4058863005744 | 3000316207 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069073 | 4058863005744 | 3000303595 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069073 | 4058863005744 | 3000318845 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069073 | 4058863005744 | 3000317551 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069073 | 4058863005744 | 3000318848 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069073 | 4058863005744 | 3000318849 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069073 | 4058863005744 | 3000319744 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069073 | 4058863005744 | 3000319745 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069073 | 4058863005744 | 3000319746 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069073 | 4058863005744 | 3000319747 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069073 | 4058863005744 | 3000319774 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069073 | 4058863005744 | 3000324501 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069076 | 4058863078298 | 3000323402 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Spain | 701069073 | 4058863005744 | 3000325551 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|--------|-----------|---------------|------------|------------------------------------|
| Spain | 701069073 | 4058863005744 | 3000325554 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069076 | 4058863078298 | 3000327600 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Spain | 701069073 | 4058863005744 | 3000301227 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069073 | 4058863005744 | 3000298132 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069073 | 4058863005744 | 3000327125 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069073 | 4058863005744 | 3000318847 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Spain | 701069073 | 4058863005744 | 3000336493 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Sweden | 701069073 | 4058863005744 | 3000281883 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Sweden | 701069073 | 4058863005744 | 3000303594 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Sweden | 701069073 | 4058863005744 | 3000317533 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Sweden | 701069073 | 4058863005744 | 3000317549 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Sweden | 701069073 | 4058863005744 | 3000316210 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Sweden | 701069073 | 4058863005744 | 3000318849 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Sweden | 701069073 | 4058863005744 | 3000319745 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Sweden | 701069073 | 4058863005744 | 3000319747 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Sweden | 701069073 | 4058863005744 | 3000319776 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Sweden | 701069073 | 4058863005744 | 3000324502 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Sweden | 701069073 | 4058863005744 | 3000324504 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Sweden | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Sweden | 701069073 | 4058863005744 | 3000324505 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Sweden | 701069073 | 4058863005744 | 3000327125 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Sweden | 701069073 | 4058863005744 | 3000318847 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Sweden | 701069073 | 4058863005744 | 3000336494 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|-------------|-----------|---------------|------------|------------------------------------|
| Switzerland | 701069073 | 4058863005744 | 3000303593 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Switzerland | 701069073 | 4058863005744 | 3000316207 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Switzerland | 701069073 | 4058863005744 | 3000303595 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Switzerland | 701069073 | 4058863005744 | 3000298132 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Switzerland | 701069073 | 4058863005744 | 3000318848 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Switzerland | 701069073 | 4058863005744 | 3000319744 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Switzerland | 701069073 | 4058863005744 | 3000319746 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Switzerland | 701069073 | 4058863005744 | 3000319774 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Switzerland | 701069073 | 4058863005744 | 3000323079 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Switzerland | 701069076 | 4058863078298 | 3000323402 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Switzerland | 701069073 | 4058863005744 | 3000324505 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Switzerland | 701069076 | 4058863078298 | 3000327124 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Switzerland | 701069073 | 4058863005744 | 3000325553 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Switzerland | 701069076 | 4058863078298 | 3000327600 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| Switzerland | 701069073 | 4058863005744 | 3000327125 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Switzerland | 701069073 | 4058863005744 | 3000318847 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Switzerland | 701069073 | 4058863005744 | 3000336493 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Thailand | 701069073 | 4058863005744 | 3000301226 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Thailand | 701069073 | 4058863005744 | 3000317533 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Thailand | 701069073 | 4058863005744 | 3000318849 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Thailand | 701069073 | 4058863005744 | 3000316210 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Thailand | 701069073 | 4058863005744 | 3000319745 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Thailand | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|----------------------|-----------|---------------|------------|---|
| United Arab Emirates | 701069073 | 4058863005744 | 3000282261 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| United Arab Emirates | 701069073 | 4058863005744 | 3000282259 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| United Arab Emirates | 701069073 | 4058863005744 | 3000282262 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| United Arab Emirates | 701069073 | 4058863005744 | 3000282260 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| United Arab Emirates | 701069073 | 4058863005744 | 3000303593 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| United Arab Emirates | 701069073 | 4058863005744 | 3000318845 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| United Arab Emirates | 701069073 | 4058863005744 | 3000318848 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| United Arab Emirates | 701069073 | 4058863005744 | 3000319744 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| United Arab Emirates | 701069073 | 4058863005744 | 3000319746 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| United Arab Emirates | 701069073 | 4058863005744 | 3000325551 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| United Arab Emirates | 701069073 | 4058863005744 | 3000325553 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| United Kingdom | 701069073 | 4058863005744 | 3000281237 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| United Kingdom | 701069073 | 4058863005744 | 3000317533 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| United Kingdom | 701069073 | 4058863005744 | 3000317549 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| United Kingdom | 701069073 | 4058863005744 | 3000318845 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| United Kingdom | 701069073 | 4058863005744 | 3000319775 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| United Kingdom | 701069076 | 4058863078298 | 3000323402 | BE-HLS 5050 #SHLS Set Advanced 5.0 |
| United Kingdom | 701069073 | 4058863005744 | 3000325554 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000280906 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000280907 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000280908 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000280909 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000282253 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|---------------|-----------|---------------|------------|---|
| United States | 701069078 | 4058863080383 | 3000282254 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000282255 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000300688 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069077 | 4058863076355 | 3000301222 | BEQ-HLS 5050 USA #SHLS Set Advanced 5.0 |
| United States | 701069078 | 4058863080383 | 3000317531 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000314506 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000316206 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000316600 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000317532 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000318851 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000318852 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000318853 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000318850 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000319748 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000319750 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000319749 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000319751 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069077 | 4058863076355 | 3000319782 | BEQ-HLS 5050 USA #SHLS Set Advanced 5.0 |
| United States | 701069078 | 4058863080383 | 3000291127 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000319778 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000319781 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000323082 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000323084 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

| | | | | |
|---------------|-----------|---------------|------------|---|
| United States | 701069078 | 4058863080383 | 3000323085 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000324495 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000324496 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000324497 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000324499 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000325557 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000325558 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069077 | 4058863076355 | 3000327602 | BEQ-HLS 5050 USA #SHLS Set Advanced 5.0 |
| United States | 701069078 | 4058863080383 | 3000325559 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000327120 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000327121 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000328404 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000328405 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000334413 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000334414 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| United States | 701069078 | 4058863080383 | 3000336487 | BEQ-HLS 7050 USA #SHLS Set Advanced 7.0 |
| Vietnam | 701069073 | 4058863005744 | 3000301226 | BE-HLS 7050 #SHLS Set Advanced 7.0 |
| Vietnam | 701069073 | 4058863005744 | 3000325555 | BE-HLS 7050 #SHLS Set Advanced 7.0 |

Template: CP-SOP-001-T-02 V02, Effective date 2019-09-15

Annex II Further information regarding Hazardous situation, Harms and Risk Levels

This Annex II Further information regarding Hazardous situation, Harms and Risk Levels is considered a supplementary attachment to the 946519 Field Safety Notice.

| Hazardous situation | Harm | S from part III | P from above | Risk | | |
|---|--------------|--------------------|-----------------|------|-----|------|
| | | | | Low | Med | High |
| Patient is exposed to pathogenic agents | Inflammation | 3 | 3 | □ | ⊗ | □ |
| | Infection | 4 | 3 | □ | ⊗ | □ |
| | Sepsis | 4 | 3 | □ | ⊗ | □ |

Severity Definitions:

Negligible (1) Inconvenience or temporary discomfort of patient, user or third party. No medical intervention or follow-up treatment is required

Low (2) Temporary injury or disability of patients, users or third parties. No medical intervention or follow up treatment is required.

Critical (3) Temporary injury or disability of patients, users or third parties. Medical intervention or follow-up treatment is required.

Catastrophic (4) Permanent injury or disability (e.g., loss of a body part), a life-threatening situation or death of patients, users or third parties

Probability Definitions:

Improbable (1) Harm is not likely.

Remote (2) Harm occurs infrequently

Occasional (3) Harm may occur occasionally / intermittent

Probable (4) Harm may occur often

Frequent (5) Harm will occur repeatedly